

Product Catalog

Capillary Lancets | Surgical Blades and Scalpels | Insulin Pen Needles

Mission

To improve health care with safe products at reduced cost

Vision

To build a mature and stable company as a respected leader in medical manufacturing

Core Values

Focus Ingenuity Innovation Integrity



SteriLance Medical (USA) Inc.
337 Curie Drive, Alpharetta, GA 30005-2264
Tel: +1(770)740-1401 Fax: +1(678)254-1771

SteriLance Medical (Suzhou) Inc.
No.168 PuTuoShan Road, New District, Suzhou 215153, China
Tel: +86 (512) 65799308 Fax: +86 (512) 67217663

SteriLance, SteriHeel, Index, Asurity and Easydrip
are the registered or pending trademarks of SteriLance Medical.
International Patents Pending
www.sterilance.com
V1.0

Since 1993

Established in 1993, SteriLance Medical is dedicated to developing, manufacturing and marketing of medical consumables. SteriLance has been making continuous efforts to provide innovative safety-engineered medical products to hospital and home care market, and now it is one of the leading manufacturers of capillary blood lancing devices and disposable scalpels worldwide.



STERILANCE™ Press Plus & Press2 Plus Pressure-Activated Disposable Safety Lancets

New double-spring design

Adequate blood volume, minimal pain



Color	Press Plus Reorde#	Color	Press2 Plus Reorde#	Gauge	Penetration Depth
	05-233015		05-243015	30G	1.5mm
	05-232818		05-242818	28G	1.8mm
	05-232618		05-242618	26G	1.8mm
	05-232122		05-242122	21G	2.2mm
	05-232118		05-242118	21G	1.8mm
	05-231720B		05-241720B	17G Blade	2.0mm

The larger the gauge number of the product, the thinner the steel lancet;
The smaller the number, the thicker the steel lancet.

Standard packing quantity: 100 pcs / box and 1,000 pcs / case.

Please contact SteriLance Medical for customized packing inquiries.

Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Refer to the instructions for contraindications or precautions.



STERILANCE™ ImPress Pro & ImPress Pressure-Activated Disposable Safety Lancets



Recommended for blood sugar monitoring devices

Tactile Contact

Stable and accurate positioning to collect an adequate volume of blood

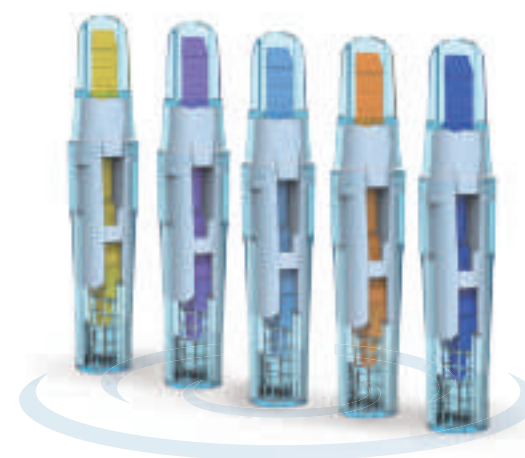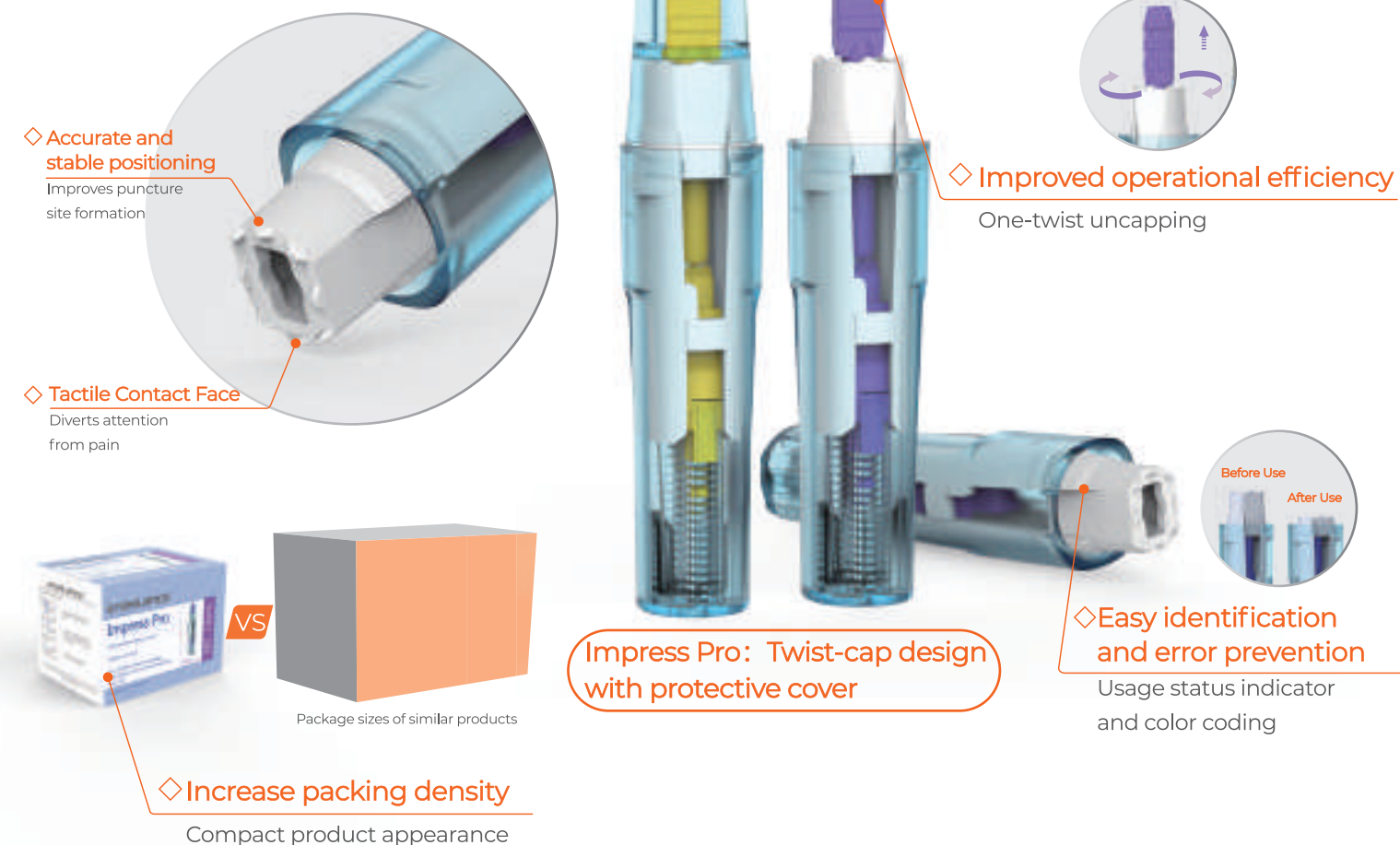

Color	ImPress Pro Reorde#	Color	ImPress Reorde#	Gauge	Penetration Depth
	05-193015		05-173015	30G	1.5mm
	05-192818		05-172818	28G	1.8mm
	05-192618		05-172618	26G	1.8mm
	05-192118		05-172118	21G	1.8mm
	05-192122		05-172122	21G	2.2mm

The larger the gauge number of the product, the thinner the steel lancet;
The smaller the number, the thicker the steel lancet.


Standard packing quantity: 100 pcs / box and 2,000 pcs / case.

Please contact SteriLance Medical for customized packing inquiries.

Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Refer to the instructions for contraindications or precautions.

STERILANCE™ Lite3

Side button-activated Disposable Safety Lancet

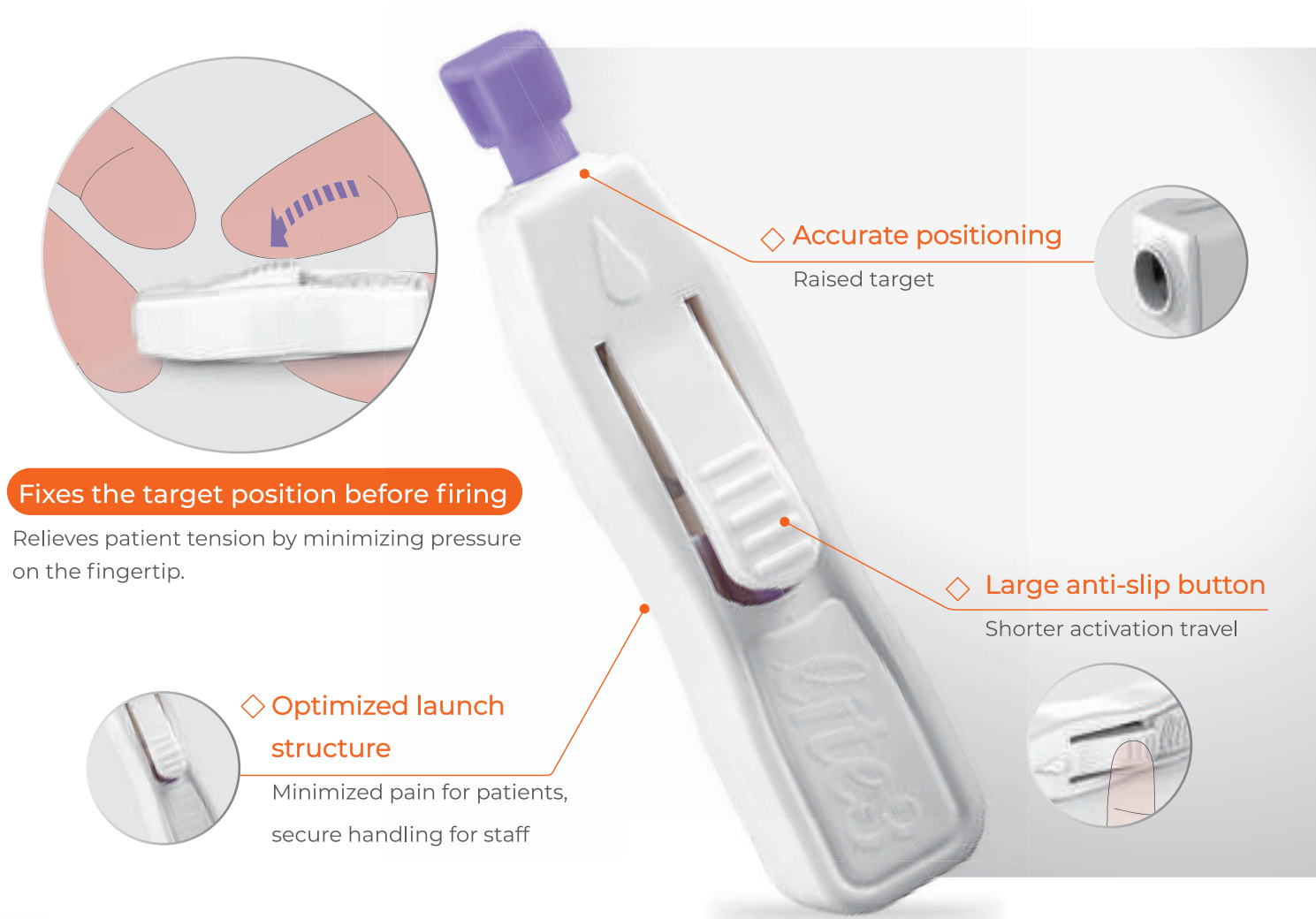

Color	Reorder#	Gauge	Penetration Depth
	05-113018	30G	1.8mm
	05-112818	28G	1.8mm
	05-112618	26G	1.8mm
	05-112118	21G	2.4mm
	05-111720B	17G Blade	2.0mm

The larger the gauge number of the product, the thinner the steel lancet;
The smaller the number, the thicker the steel lancet.

Standard packing quantity: 100 pcs / box and 1,000 pcs / case.
Please contact SteriLance Medical for customized packing inquiries.
Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Refer to the instructions for contraindications or precautions.

STERILANCE™ Flex3

Depth-adjustable, Pressure-activated lancing Disposable Safety Lancet



Color	Reorder#	Gauge	Penetration Depth
	05-1528A	28G	1.3/1.8/2.3mm
	05-1526A	26G	1.3/1.8/2.3mm
	05-1523A	23G	1.3/1.8/2.3mm


The larger the specification number of the product, the thinner the steel lancet;
The smaller the number, the thicker the steel lancet.

Standard packing quantity: 100 pcs / box and 1,000 pcs / case.
Please contact SteriLance Medical for customized packing inquiries.
Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Refer to the instructions for contraindications or precautions.

STERILANCE™ Elite


Rear button activated Disposable Safety Lancet

Gamma-sterilized, assuring sterility for the patient

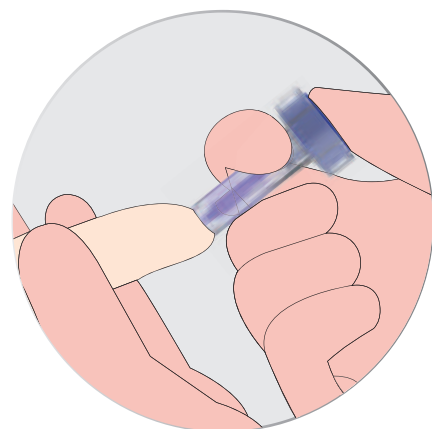

Safer T-shape design

Best fit for nurses needing gentle skin contact
Large anti-slip rear button providing a stable and comfortable pressing surface.

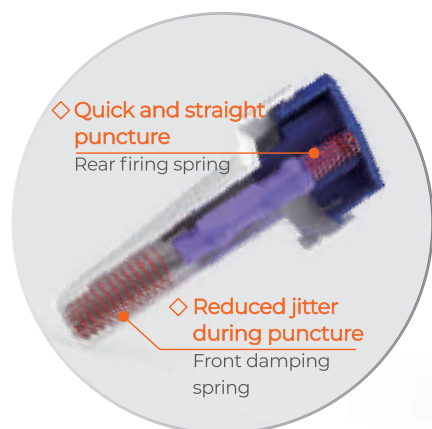

Comfortable double-spring action



Color	Reorder#	Gauge	Penetration Depth
	05-213015	30G	1.5mm
	05-212818	28G	1.8mm
	05-212618	26G	1.8mm
	05-212118	21G	1.8mm
	05-211720B	17G Blade	2.0mm

The larger the gauge number of the product, the thinner the steel lancet;
The smaller the number, the thicker the steel lancet.

Standard packing quantity: 100 pcs / box and 1,000 pcs / case.


Please contact SteriLance Medical for customized packing inquiries.

Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Refer to the instructions for contraindications or precautions.

STERILANCE™ Elite Pro

Depth-adjustable, Rear button activated Disposable Safety Lancet

Single SKU helps simplify inventory management

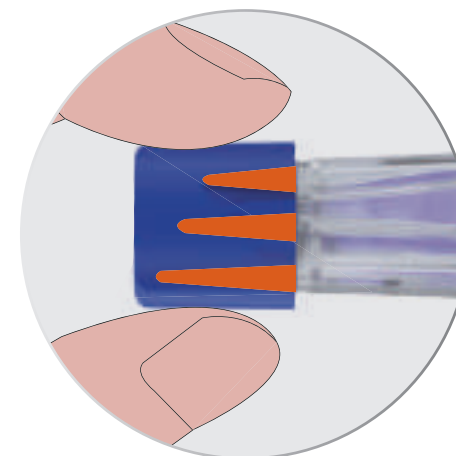

One Lancet offers three depth settings (1.3, 1.8, and 2.3mm) for various patient conditions

Simply select the depth that provides the greatest comfort—while at the patient's bedside—and easily obtain an adequate blood sample for single dose tests.

Trigger



Color	Reorder#	Gauge	Penetration Depth
	05-2928	28G	1.3/1.8/2.3mm
	05-2926	26G	1.3/1.8/2.3mm
	05-2923	23G	1.3/1.8/2.3mm
	05-2921	21G	1.3/1.8/2.3mm

The larger the specification number of the product, the thinner the steel lancet;
The smaller the number, the thicker the steel lancet.

The standard packing quantity: 100 pcs / box and 1,000 pcs / case.

Please contact SteriLance Medical for customized packing inquiries.

Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Refer to the instructions for contraindications or precautions.



Newborn disease screening

Specialized Heel Incision Safety Lancets




SteriHeel™ Plus

Heel Incision Safety Lancets

◇ One-handed blood drawing to increase efficiency

The one-piece safety lock can be released with one hand
Position the heel with one hand and push the button to fire

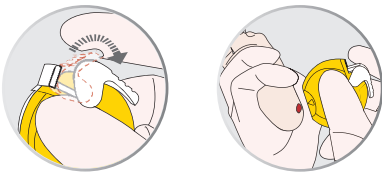

◇ Micro-blade with patented gentle launch mechanism and smooth blade path

Precise depth control and adequate blood volume without massaging reduces hemolysis, improves test accuracy and speeds the small incision healing




◇ Specialized product for newborns with assured quality and safety

Individually blister packaged and gamma radiation-sterilized
Blade retracts after use, preventing cross-contamination




Neoheel Heel Incision Safety Lancets

Ordinary safety lancets

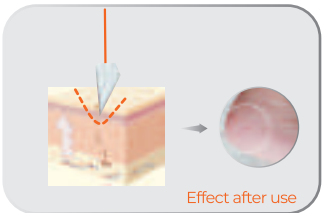

Shallow incision makes healing faster

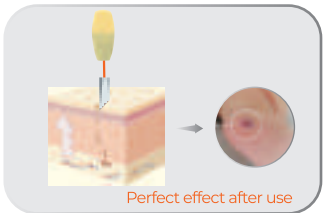

Deep incision makes healing slower

VS

Color	Reorder#	Incision Depth	Incision Length
	30-106514	0.65 mm	1.4mm
	30-101025	1.0 mm	2.5mm
	30-108575	0.85 mm	1.75mm
	30-102030	2.0 mm	3.0mm

Standard packing quantity: 50 pcs / box and 200 pcs / case.
Please contact SteriLance Medical for customized packing inquiries.
Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Refer to the instructions for contraindications or precautions.

Improving diabetes management
Helping diabetics live a better life

Creating a "pain-free" blood sampling experience
Enhancing blood sugar monitoring

The perfect partner for home-use blood sugar test strips



STERILANCE™ LDE4

Adjustable Lancing Device with Ejector

◇ Patented SSC (Soft Stop Control) technology absorbs shock and reduces pain



Soft-Stop Control

At the moment of puncture, the needle engages the SSC structure which absorbs the impact of the needle and stabilizes the puncture, through elastic deformation and energy conversion.



◇ 9 adjustable depth levels

Easy to read settings



◇ Optional transparent AST cap



Patented inner cap cleaning technology

Traditional caps consist of inseparable inner and outer layers. During long-term lancet use, blood residue may collect on the inner cap surface, risking bacterial growth. The LDE4 features a cleaning setting (↕ Setting). At this setting, the outer cap can be removed, to enable easy cleaning of the inner cap.



Rotate to the "↕" setting



Remove the outer cap



Wipe and clean

Reorder#	Description	Color	Package
02-2201	9 Adjustable depths	Gray/White	Individually packaged 100 pcs/case



Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Refer to the instructions for contraindications or precautions.

- ◇ Comfort grip for needle placement
- ◇ International design benchmarked
- ◇ Well-established reputation for quality



Pro cap design





















Safe disposal of fine needles



Softsure

protection against fine
needle deformation during loading



Color	Soft Reorder#	Color	Softsure Reorder#	Color	Soft Pro Reorder#	Color	Softsure Pro Reorder#	Gauge
	01-0132		01-1732		01- 0232		01-1832	32G
	01-0130		01-1730		01-0230		01-1830	30G
	01-0128		01-1728		01- 0228		01-1828	28G
	01-0126		01-1726		01- 0226		01-1826	26G
	01-0121		01-1721		01-0221		01-1821	21G

The larger the gauge number of the product, the thinner the steel lancet;
The smaller the number, the thicker the steel lancet.


Standard packing quantity: 100 pcs / box and 10,000 pcs / case.

Please contact SteriLance Medical for customized packing inquiries.

Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Refer to the instructions for contraindications or precautions.

Alcohol Prep Pads


- Isopropyl Alcohol, 70% by volume.
- For preparation of the skin prior to injection or lancing.
- For external use only.

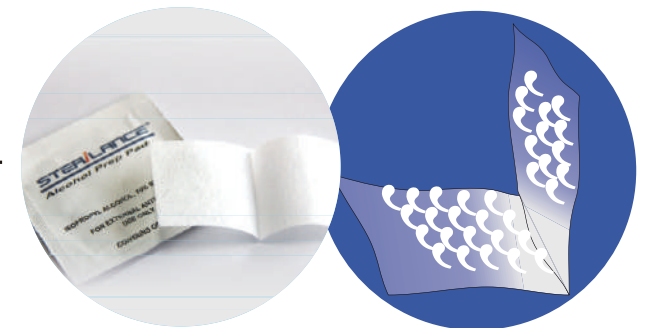

Reorder #	Description	Specification	Package
06-023060	Alcohol pre pads	30*60mm	1pc/sealed foil wrap 100pcs/box 5,000pcs/case

Skin Graft

Designed for dedicated skin graft handles.Individually packaged and gamma sterilized, with a 5 year shelf life.

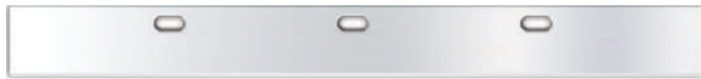

Reorder#	Description	Gauge	Package
13-01158	Carbon Steel	158mm	1 pc/sealed foil wrap 50 pcs/box 1,000 pcs/case
19-01158	Stainless Steel		

Stitch Cutter

A rounded bull-nose tip design for stitch removal while minimizing pressure and damage to the skin.Individually packaged and gamma sterilized, with a 5 year shelf life.



Reorder#	Description	Gauge	Package
14-01XX	Carbon Steel	No. 65#,110#	1 pc/sealed foil wrap 100 pcs/box 5,000 pcs/case
15-01XX	Stainless Steel		

Gouge Blades

Suitable for podiatry and foot treatment, such as corns, calluses and hard skin removal.Precision blade profiles for specific areas and operations with easy interchangeability.Individually packaged to meet hygiene standards, with a half year service life for carbon steel and one year for stainless steel.




Reorder#	Description	Gauge	Package
24-01XX	Carbon Steel	No. 1#,2#,3#,4#,5#, 6#,8#,10#,12#,15#	1 pc/sealed foil wrap 100 pcs/box 5,000 pcs/case
24-02XX	Stainless Steel		

Notes: XX represents the blade specification.
Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Kindly refer to the manual for details about contraindications and precautions


Disposable Scalpel

◇ Patented safety lock

Efficient one-handed sliding sheath action

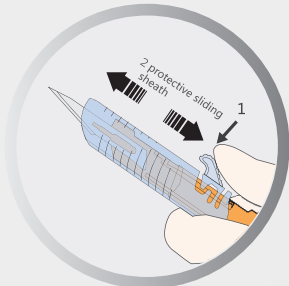

◇ Anti-reuse lock quickly activated after use

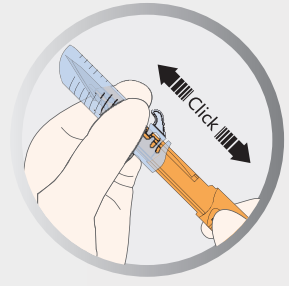

◇ #3 and #4 handles
common blade specifications



◇ Suitable for both right- and left-handed use
Top button design

◇ On-site measurement
Scales in cm

◇ Integrated blade and sheath design reduces the risk of cross-infection

Protective sheath slides to expose/cover the blade, reducing potential accidental stab injuries to healthcare professionals during attachment and removal of blades as well as passing of scalpels

◇ Fast operation and high visibility
common blades pre-installed

◇ Individually packaged
Blister packaged , Gamma sterilized
5 year shelf life



Reorder#	Description	Gauge	Package
22-03XX	Index Asurity	No . 9#-17#	1 pc/blister packed 10 pcs/box 500 pcs/case
22-04XX		No . 18#-36#	

Notes: XX represents the blade specification.
Please read the product instructions carefully or purchase and use under the guidance of a medical professional. Kindly refer to the manual for details about contraindications and precautions



Disposable Scalpel (Standard I)

◆ **Efficient use**
Blade pre-installed with a full range of blades available

◆ **#3 and #4 handles**
common blade specifications

◆ **Anti-slip grip**

◆ **Facilitating on-site measurement of incision size and lesion diameter**
Scales in cm

◆ **Blade and cover design reduces the risk of cross-infection**
Removal / exposure of the protective cover exposes/covers the blade, reducing potential accidental stab injuries to healthcare professionals during attachment and removal of blades as well as passing of scalpels

◆ **Individually packaged**
Blister packaged, Gamma sterilized
5 year shelf life

Reorder#	Description	Gauge	Package
07-01XX	Standard I	No . 9#-17#	1 pc/blister packed 10 pcs/box 500 pcs/case
07-02XX		No . 18#-26#	

Notes: XX represents the blade specification.
Please read the product instructions carefully or purchase and use under the guidance of a medical professional.
Kindly refer to the manual for details about contraindications and precautions



Disposable Scalpel (MINI)

◆ Blade pre-installed and ready to use upon opening

◆ Small and light-weight, portable

◆ Customized packaging available to suit customer needs

Reorder#	Description	Gauge	Package
16-01XX	MINI	No . 9#-17#	40 pcs/box 2,000 pcs/case

Notes: XX represents the blade specification.
Please read the product instructions carefully or purchase and use under the guidance of a medical professional.
Kindly refer to the manual for details about contraindications and precautions



Disposable Surgical Blades

◆ Available in carbon steel and stainless steel, meeting different needs

◆ No. 9-17 blades fit No. 3 scalpel handles, No. 18-36 blades fit No. 4 scalpel handles

◆ **Single-piece sealed**
Individually packaged, and sterilized by gamma ray irradiation, with a 5 year shelf life

◆ **Consistently sharp cutting edge for smooth and clean cutting**

◆ **Full range of blade specifications available with standardized handle fit.**

Reorder#	Description	Gauge	Package
10-01XX	Carbon Steel	No .9#-17# No.18#-36#	1 pc/foil seal 100 pcs/box 5,000 pcs/case
11-01XX	Stainless Steel	PM40#、 60#、 70#	

Notes: XX represents the blade specification.
Please read the product instructions carefully or purchase and use under the guidance of a medical professional.
Kindly refer to the manual for details about contraindications and precautions

Company Introduction

Sandstone Medical is a state of the art medical device manufacturer located in Suzhou Science and Technology Town. Sandstone offers project management, concept brainstorming, prototype testing and validation, and the knowledge to traverse complex and emerging regulatory compliance. Sandstone follows the strictest inspection procedures with the most advanced equipment available, including continuous improvement software and intelligent eye inline cameras to ensure the highest output.

Sandstone not only produces the products but we help OEMs and contract manufacturers get their devices to market within tight development timelines while meeting strict ISO standards. When your company needs a complete source for highly engineered products, Sandstone delivers excellence.



MDSAP

MDL

510K

CE

FMR

Dedicated Manufacturer



Innovation

Constant quality improvement, innovation and new technologies are the core of our heritage.



Dedicated

Specialized in development and manufacturing of needle consumables.

Quality Assured



Precision

High precision injection machines, full automatic assembly lines, packaging machines from German, Japan and Korea.



Professional

Product on-line automatic testing by carema system to ensure consistent quality.

easydrip™ Classic

Pen Needles



32G	31G	31G
4mm	5mm	6mm

*Pen injectors compatible with the pen needle:
NovoPen®5, Autopen® 24, Gansulin Pen, Ypsopen®, UNIPEN®, XiuLin Pen,
HumaPen®ErgoII, Byetta, Victoza Pen
* All the brands are trademarks of their respective owners.
Please read insulin pen instructions for use.






High Flow Rate

Silicone Coating

Classic bevel

Gamma Irradiation

- Sterilized by Gamma Irradiation
- Non-Toxic, Non-Pyrogenic
- Without Natural Rubber Latex
- Full range sizes, Private Label Available

Ordering information				
Reorder #	Color	Needle gauge	Needle length	Standard packaing
01-3204		32G needle	4 mm	100 pcs/box 2,000 pcs/case
01-3105		31G needle	5 mm	
01-3106		31G needle	6 mm	

1. Needle gauge refers to the size of the needle. The higher the gauge, the thinner the needle;

2. Please contact our sales person if you have special requests for an unique size.

easydrip™ Plus

Pen Needles

Innovation for easy & conformable injections

New generation pen needles provide with more comfort against the skin during the injection process.



Reliable Injection Depth



High Flow Rate




Silicone Coating



New Generation Design

More comfortable against the skin during the injection process

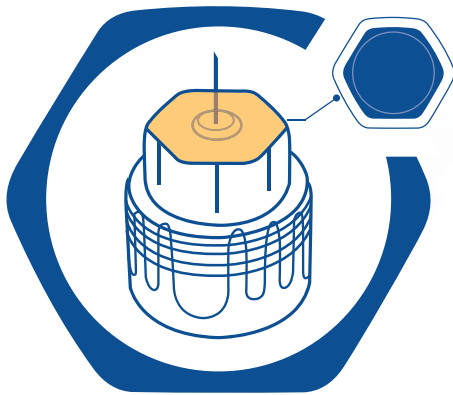


*Pen injectors compatible with the pen needle:
 NovoPen®5, Autopen®24, Gansulin Pen, Ypsopen®, UNIPEN®, XiuLin Pen,
 HumaPen®ErgoII, Byetta, Victoza Pen
 * All the brands are trademarks of their respective owners.
 Please read insulin pen instructions for use.

Contoured hub for Reliable Depth Control

The contoured hub distributes insertion forces across a wider skin surface with the intent to minimize the clinical impact (IM injection) of variable patient injection forces on injection depot depth.



Intramuscular (IM) injection can increase insulin absorption, causing hypoglycemia.

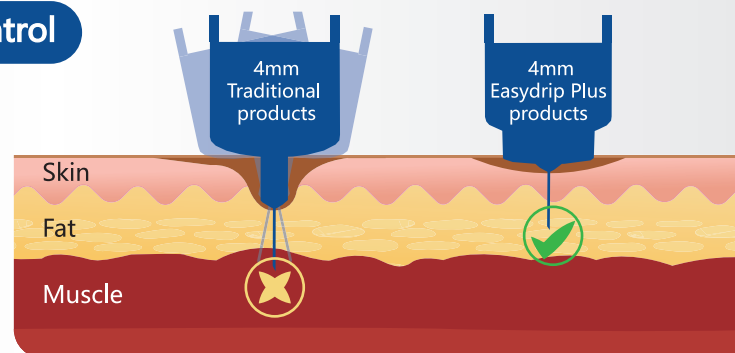


Ultra-strong Design to resist bending or breaking



An improved ability to hold the pen securely against the skin without wobbling during an injection.

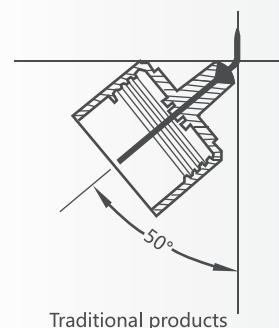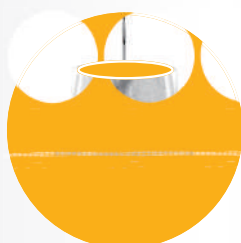


Traditional products

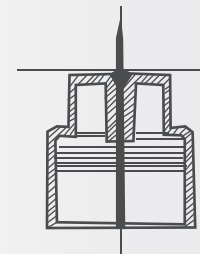

Easydrip Plus products

Ergonomic design helping diabetes use PN's more easily

Wider outer cover with ergonomically designed grooves and a colored, larger, more prominent inner shield makes it easier to remove the inner shield.

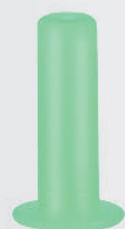

Wider inner cap



Wider outer cap



Ordering information

Reorder #	Color	Needle gauge	Needle length	Standard packaing
02-3204		32G needle	4 mm	100 pcs/box 2,000 pcs/case
02-3105		31G needle	5 mm	
02-3106		31G needle	6 mm	

1. Needle gauge refers to the size of the needle. The higher the gauge, the thinner the needle;
2. Please contact our sales person if you have special requests for an unique size.

Sandstone

Keep safety insulin injections short and simple




easydrip™ ReadyGard™ Duo
Safety Pen Needles
With Dual Automatic Protective Shields

easydrip™ ReadyGard™ Duo

Safety Pen Needles HS Model

With Dual Automatic Protective Shields

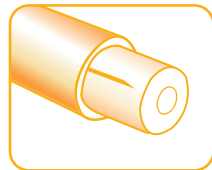


Reliable Injection Depth



High Flow Rate



Silicone Coating

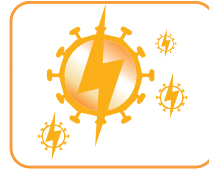

Gamma Irradiation

Safety-Engineered Pen Needle with Passive Dual-Ended Protection



Chinese patents granted and US and European patents pending

*Pen injectors compatible with Safety Pen Needle:
NovoPen®5, Autopen®24, Gansulin Pen, Ypsopen®, UNIPEN®, XiuLin Pen,
HumaPen®ErgoII, Byetta, Victoza Pen
* All the brands are trademarks of their respective owners.
Please read insulin pen instructions for use.

Diabetes nurses at risk of blood-borne pathogens

In a survey performed amongst 634 European diabetes nurses in 2012, nearly a third reporting suffering a NSI in the context of giving injections to patients with diabetes in the hospital setting, and approximately 24% of US nurses have suffered a NSI from giving diabetic injections.



NSI occur with both ends, and both can be sources of potential infection

Pen injection devices aspirate human cells back into the cartridge during use. These potentially infectious cells can then be returned to the needle and accidentally transmitted to another person via a NSI. Both ends of the pen needle can do this. European data has shown that up to 10% of NSI with pen needles occur with the cartridge end, not the patient end. Hence protection must be provided against NSI from BOTH ends of the pen needle.

Easydrip ReadyGard Duo with Passive, Dual-Ended Protection

Before Use


On the patient end, The needle itself is guarded within a plastic sheath allowing it to be undetectable until it is ready to be used.

On the cartridge end, the sterile needle is exposed to connect the injection devices



On the patient end, a red indicator confirms the safety mechanism has been activated.

On the cartridge end, protection is confirmed when the red shield covers the needle upon removal from the pen.

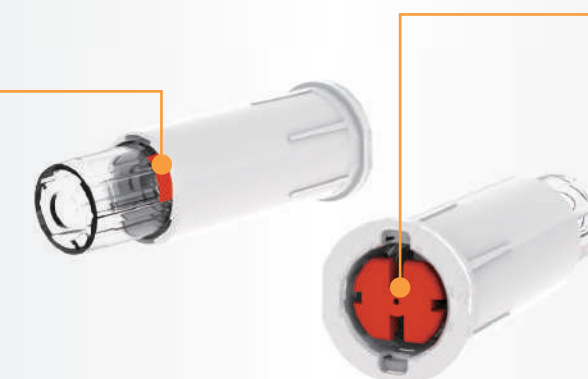


After Use

Safety devices such as a pen needle with this dual protection have been shown to reduce the incidence of NSI dramatically, and accounted for <2% of NSI in a survey.

Less steps improving efficiency

No inner cap equipped and no need to cover back outer cap before discarding the used needles.

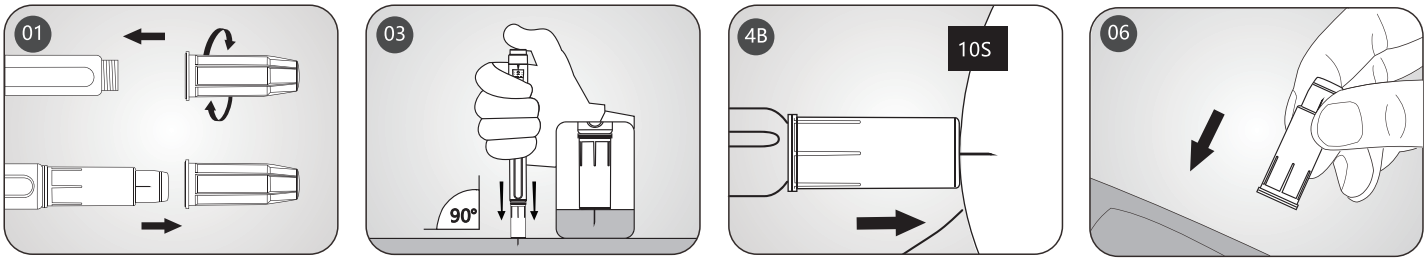


Comfort injection

The patented structure allows for a reliable and straight needle in and out path, hence to provide the patient with a more comfortable and smooth injection experience.

Controllable depth

A wider skin surface, not only improve patient comfort due to less pressure during injection, also with the intent to minimize the clinical impact (IM injection) of variable patient injection forces on injection depot depth.

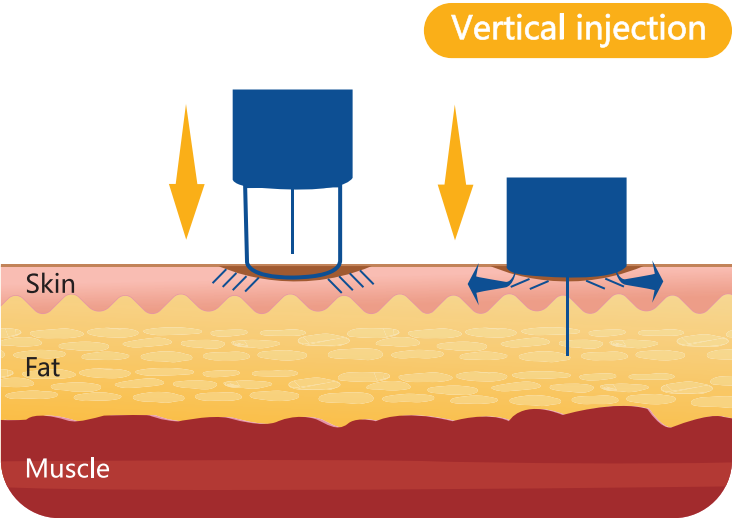

Prevent the injection-related anxiety

The Readygard Duo Needle has a front shield that covers the needle to help prevent the injection-related anxiety before and after the injection. It has been noted by pharmacists and prescribers alike that the less a person sees a needle, the less stress and anxiety they will have.



Ordering information

Reorder #	Color	Needle gauge	Needle length	Standard packaing
05-3105		31G needle	5 mm	100 pcs/box 800 pcs/case
05-3008		30G needle	8 mm	
05-3005		30G needle	5 mm	

1. Needle gauge refers to the size of the needle. The higher the gauge, the thinner the needle;
2. Please contact our sales person if you have special requests for an unique size.